

# STP-320 Category B Shipping System, Insulated

## **Description:**

The STP-320 is an Insulated Category B Shipping System designed for the shipment of temperature controlled shipments using phase change materials such as the SAF-T-TEMP<sup>™</sup> line as well as Dry ice as an energy source.

### STP-320 Category B, Insulated Shipping System

Quantity - 4/Case

#### System Components:

- Outer Box OD 15 3/8 x 13 5/8 x 12 3/4"
- Insulated Chest ID 11 5/8 x 10 x 8"
- STP-111 Inner Box ID 6 1/2 x 6 1/2 x 8"
- **STP-710** Disposable 2-part Secondary Pressure Vessel, Medium (2)
- STP-152 250 mL Absorbent Strip (2)
- STP-610 12 x 12" Saf-T-Rap® Bubble Wrap (2)
- STP-804 Dry ice Marking
- STP-803 Class 9 Hazard Label
- STP-818 UN 3373 Marking



Maximum Capacity - Two 500 mL blood bags or thirty-two 10 mL vials.

Ideal for long distance frozen shipments.

Designed for use with Saf-T-Temp<sup>™</sup> products or Dry ice for thermal control.


#### **Refurbishments and Optional Components**

STP-321 Outer Box - 10/Case STP-111 Inner Box - 12/Case

Specifications subject to change without notice



## STP-320 Category B Shipping System, Insulated

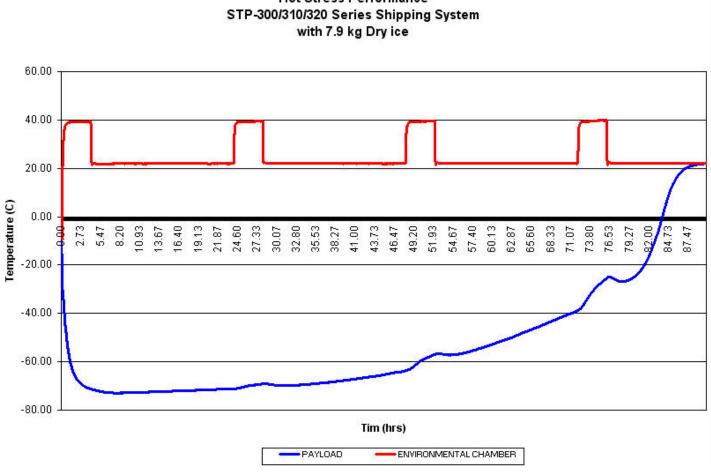

# Hot Stress Performance

#### **Testing Synopsis**

Product performance testing under worst-case conditions was done on the STP-320 Category B Insulated Shipping System. Seven and nine tenths kilograms (7.9 kg) of Dry ice was used as refrigerant.

The minimum payload of five (5) 10 mL vials filled with a non-freezing liquid mixture were placed in the STP-710 Disposable Secondary Pressure System (included). The samples were preconditioned in a conventional freezer set to maximum between -25 to -30 °C for a minimum of 24 hours prior to testing.

The average total test time where the minimum payload samples were at or below 0 °C was 83.9 hours.

The thermal profiles displayed in this data sheet are applicable to the STP-300 Series shippers, which include the STP-300; STP-310; and STP-320.